

# RELAX

# SNUG FIT BACK



Rear Back Angle



Inner Foam and Incontinent Cover

The **Relax Snug Fit** by Power Plus Mobility is designed to meet a variety of positioning needs. This is an aluminum light weight formed back with hardware that allows for multiple adjustments. The Relax Snug Fit with its open hip laterals is suited for users needing complex postural support and spinal alignment. Provides 7" of mid Thoracic support.

### Quick Reference Guide

#### RELAX SERIES - SNUG FIT BACK

**REGULAR SIZES ADP: [SEBWM0027](#)**

Width = 16" - 18"

Part Number: Refer below

Height = 18"

**CUSTOM ADP: WIDTH [SEBWM0035](#), HEIGHT [SEBWM0040](#)**

Width = 15" or <, 19" or >

Part Number: Refer below

Height = 17" or <, 19" or >

"<" means less and ">" means more


#### PART NUMBERS

Part Number: MDU(W X H)

Fill in width (W) and height (H) to get part number.


E.g. For a W = 15" and H = 17" Back; Part number would be MDU1517

# RELAX

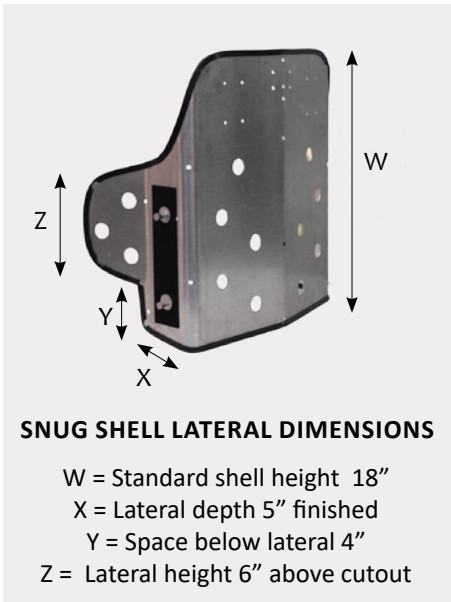

**SNUG SHELL LATERAL DIMENSIONS**

W = Standard shell height 18"  
 X = Lateral depth 5" finished  
 Y = Space below lateral 4"  
 Z = Lateral height 6" above cutout

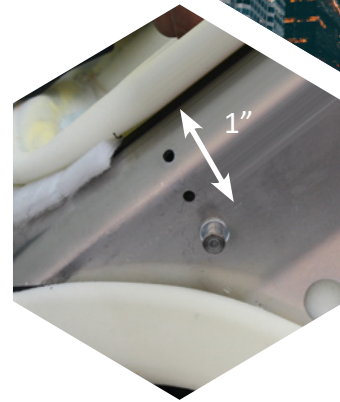


**BACK PIN POSITIONS ADJUSTABLE HOLES**

Ability to reposition pin hardware back by 1" on the shell. Optional 1" longer pins to mount a smaller back on a wider chair

**OVERLAYS**

\* Foam Pad      \* Visco      \* Gel

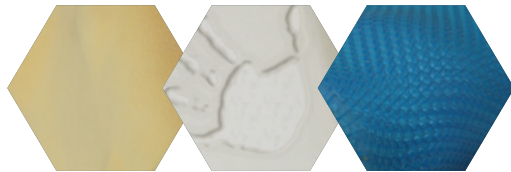

ADP: SEBND0015      SEBND0015      SEBND0020

\* Up-charge items

**KEY FEATURES**

- Provides complex posterial support with open hip laterals
- Comfortable High Resilience Foam
- 3 lbs. Aluminum, Light Weight Shell
- Pre Punched holes to mount headrest and laterals
- Adjustable/Removable Lumbar Support
- Easy to Install and Remove
- Height Adjustable



**UNIVERSAL HARDWARE**

Accommodates Pelvic rotation and will allow for full back support with 25 degrees of posterior and 20 degrees of anterior adjustability.

SPECIFICATIONS	RELAX SERIES - SNUG FIT BACK
ADP	SEBWM0027, CUSTOM W SEBWM0035, H SEBWM0040
STANDARD WIDTH RANGE	16" - 18"
CUSTOM WIDTH RANGE	17" or less and 19" or more
STANDARD HEIGHT RANGE	18"
CUSTOM HEIGHT RANGE	17" or less and 19" or more
AVAILABLE HARDWARE STYLES	One style which is universal to all our backs
ANGLE ADJUSTABLE	25° anterior adjustability, 20° of posterior adjustability
WIDTH ADJUSTABLE	2" of adjustability through the pins on the back
FITS TUBING SIZE	7/8" tubing or 1"
HEADREST MOUNT	Backs pre-drilled for headrest supports and lateral supports

SNUG FIT BACK

RELAX

BACKS