

EXTREME


TILT WHEELCHAIR



ADP CODE: WM0001041

AADL CODE: W305

The EXTREME tilt wheelchair has a light-weight frame and is depth adjustable. The tilt system offers users a change in orientation to redistribute pressure from one area to another. The front seat to floor height on the chair also remains constant during a 0 to 25 degree tilt.

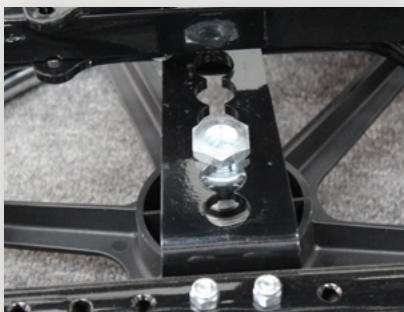

Seat To Floor Height
Adjustability



Zero Maintenance Soft Set Cylinders



Offset Front Hangers
(Optional)



Dual Cylinders Allows For Consistent Tilt For All Weight Capacities



Standard Angle Adjustable Stroller Handles



Angle Adjustable Front Journals

CUSTOMIZE ANYTHING

Power Plus Mobility's advanced machinery specializes in manufacturing custom parts to customers' specifications from either a drawing or a sample. The parts we manufacture can be simple or more complex and can be manufactured in small or large quantities. Contact your local sales representative with your specifications.

EXTREME

 **POWER PLUS MOBILITY**
OUR PRODUCTS WILL MOVE YOU

STANDARD FEATURES

- Light Weight Frame
- Depth Adjustable Frame
- ABS Seat Pan
- Dual Gas Struts for Easy Tilt (Up To 25°)
- New Recessed Lower Frame For Foot Propulsion
- Dual Gas Cylinders for Smooth & Easy Tilting
- Relax Seating System
- New Gas Cylinder Release Mechanism
- Black Frame Colour
- T-Style Arms
- Adjustable Stroller Handle
- Swing In & Out Footrest
- Adjustable Caster Journals
- New Recessed Axles
- Maintains Front STF Height In Tilt
- Overall Wheelchair Width: 29" (18" Seat Width)

SPECIFICATIONS	EXTREME
SEAT WIDTH:	14" TO 36"
SEAT DEPTH:	(ADJUSTABLE): 16" TO 24"
SEAT HEIGHT:	(ADJUSTABLE): 12.25" TO 18.25"
Tilt ANGLE:	0° TO 25°
BACK ANGLE:	(ADJUSTABLE): 0°, 5°, 10° OR 15°
ADJUSTABLE STROLLER HANDLE:	0° TO 180°
WEIGHT CAPACITY:	UP TO 600 lbs.
WEIGHT OF CHAIR:	60 lbs. (APPROXIMATELY)
	DYNAMIC TILT ADJUSTMENTS 0° TO 25°
	TWO-WAY LOCKING GAS STRUT
	SINGLE RELEASE, DUAL CYLINDERS
	STANDARD LOCKOUT COLLARS

EXTREME
TILT WHEELCHAIR



Learn more about this product by scanning the QR Code