

# FULL FORM FIT BACK



RELAX  
SERIES



The Relax Form Fit Back with its *Full* laterals is suited for users needing complex postural support and spinal alignment. This back is suited for complex seating and positioning needs. Cubes are removable to custom fit for spinal anomalies.

### Quick Reference Guide

RELAX SERIES - FULL FORM BACK		
<b>STANDARD</b> ADP: <b>SEBWM0026</b>	<b>CUSTOM WIDTH</b> ADP: <b>SEBWM0026 &amp; SEB-WM0035</b>	<b>CUSTOM HEIGHT</b> ADP: <b>SEBWM0026 &amp; SEBWM0040</b>
Height: 18" Width: 16" & 18"	Height: 18" Width: All except 16" & 18"	Height: All except 18" Width: 16" & 18"
Available Heights: 16", 18" & 20"		

**CUSTOM**  
ADP: **SEBWM0026** W= **SEBWM0035** & H= **SEBWM0040**

Part Number: Refer below

### PART NUMBERS

SHUD (W X H) **Full Laterals**

Fill in width (W) and height (H) to get part number.


E.g. For a W = 15" and H = 17" Back; Part number would be SHUD1517

### Poly Foam Cubes Part Numbers:

QUB6x6	Poly Foam Cubes 6x6
QUB12x12	Poly Foam Cubes 12x12


**OVERLAYS**

\* Foam Pad      \* Visco      \* Gel

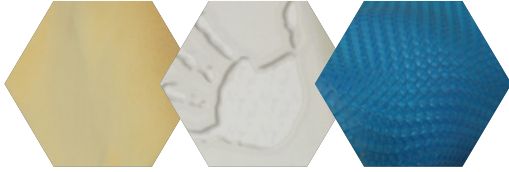

ADP:    SEBND0015    SEBND0015    SEBND0020

\* Up-charge items

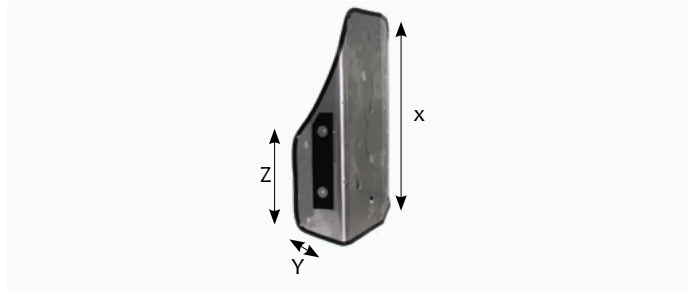

X = Standard shell height 18"  
 Y = Lateral depth 5" finished  
 Z = Lateral height 10"

**CUSTOMIZE ANYTHING**


Power Plus Mobility's advanced machinery specializes in manufacturing custom parts to customers' specifications from either a drawing or a sample. The parts we manufacture can be simple or more complex and can be manufactured in small or large quantities. Contact your local sales representative with your specifications.

**KEY FEATURES**

- Full fit back lateral with removable cubes for custom fit
- Comfortable High Resilience Foam
- 3 lbs. Aluminum, Light Weight Shell
- Pre Punched holes to mount headrest and laterals
- Easy to Install and Remove
- Height Adjustable
- TEK4 Waterproof material at the bottom of back for easy clean up



Inner Foam and Incontinent Cover

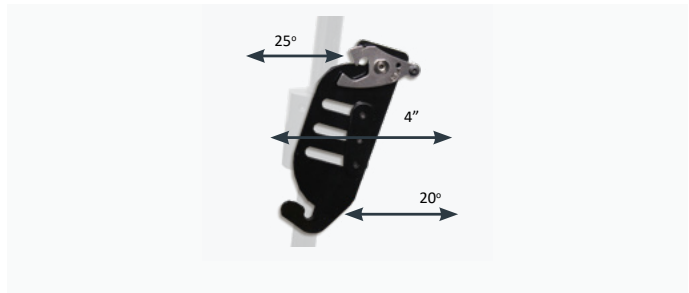

Accommodates Pelvic rotation and will allow for full back support with 25° of posterior and 20° of anterior adjustability.

SPECIFICATIONS	RELAX SERIES - FORM FIT BACK
ADP	<b>SEBWM0026 W= SEBWM0035 &amp; H= SEBWM0040</b>
STANDARD WIDTH RANGE	16" - 18"
CUSTOM WIDTH RANGE	15" or less and 19" or more
STANDARD HEIGHT RANGE	18"
CUSTOM HEIGHT RANGE	17" or less and 19" or more
AVAILABLE HARDWARE STYLES	One style which is universal to all our backs
ANGLE ADJUSTABLE	25° anterior adjustability, 20° of posterior adjustability
WIDTH ADJUSTABLE	2" of adjustability through the pins on the back
FITS TUBING SIZE	7/8" tubing or 1"
HEADREST MOUNT	Backs pre-drilled for headrest supports and lateral supports
FINISHED HEIGHT	Back Shell Size +2"